

AUTOMATIC PYCNOMETERS

THE DENSITY ANALYSIS OF POWDERS,
FOAMS AND BULK SOLIDS




Quantachrome
INSTRUMENTS

PENTAPYC 5200e

ULTRAPYC 1200e

Sample Cells	Nominal Volume	Internal Diameter	Internal Depth
Standard sample cells	135 cm ³	49 mm	75 mm
	50 cm ³	40 mm	39 mm
	10 cm ³	24 mm	23 mm
Optional Microcell	4.25 cm ³	15 mm	24 mm
Optional Mesocell	1.75 cm ³	13 mm	13 mm
Optional Nanocell	0.25 cm ³	8 mm	6.5 mm

All cells are supplied in stainless steel for durability.(aluminum available)
Non-elutriating cells are available for fine powders, and spools are available for films and fibers.

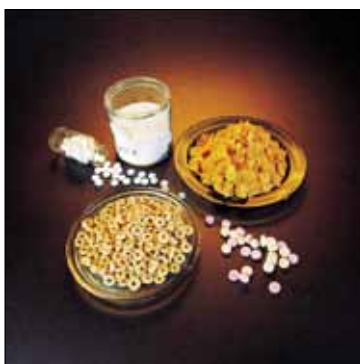

Foods and Pharmaceuticals



Catalysts and Ceramics



NIST traceable calibration spheres.

GAS PYCNOMETRY

"Pycnometry" is derived from the Greek word *puknos*, which has long been identified with volume measurements. The pycnometers from Quantachrome are specifically designed to measure the true volume of solid materials by employing Archimedes' principle of fluid displacement and gas expansion (Boyle's Law).

Ideally, a gas is used as the displacing fluid since it penetrates the finest pores assuring maximum accuracy. For this reason helium is recommended, since its small atomic dimension enables entry into crevices and pores approaching 0.2 nm. Its behavior as an ideal gas is also desirable. Other gases such as nitrogen can be used, often with no measurable difference.

Applications

Quantachrome pycnometers are used for research, development and quality control in such diverse industries as carbon black, catalysts, cement, ceramics, charcoals, cosmetics, desiccants, fertilizers, fibers, fillers, insulating and structural foams, powdered foods, ion exchange resins, minerals such as alumina, silica, titania and others, nuclear fuels, petrochemicals, pharmaceuticals and powdered metals. Gas pycnometers can even determine the percentage of solids in a slurry, or density of oil.

The **PENTAPYC 5200e** and the **ULTRAPYC 1200e** from Quantachrome are the ultimate instruments for measuring the true volume and density of powders, foams and bulk solids. A wide range of sample cell sizes are easily interchanged to accommodate many different samples. Calibration, sample conditioning, operation of valves and calculation of results are completely automatic. Samples are quickly and automatically analyzed as many times as necessary to achieve the user desired % deviation from mean at the specified number of runs. Should the deviation setting be too narrow, the analysis terminates upon reaching an operator specified maximum number of runs. The results are printed automatically, freeing the operator to perform other laboratory tasks.


The **PENTAPYC 5200e** permits up to five samples to be loaded and then purged simultaneously. Each of the five samples is then automatically analyzed in sequence with no operator involvement. Adjustable flow rates for each cell provide maximum flexibility of operation.

For laboratories with lower analytical throughput needs, the single sample station **ULTRAPYC 1200e** provides the same ease of operation and high performance as the **PENTAPYC 5200e**.

- For measurements that require a fixed, known temperature, a temperature control option is available for both pycnometers.
- For measurements in a hazardous environment, such as a nuclear "hot cell", a **remote** operation **ULTRAPYC 1200e** is available.
- For measurement of foams, specialized versions called the **ULTRAFOAM** and **PENTAFOAM** are available for automatic measurement of open and closed cell content according to ASTM D-6226 with automatic functions for the analysis of cell compressibility and cell fracture.
- For measurements of extremely small amounts of material, the **MICRO-ULTRAPYC 1200e** provides three interchangeable cell sizes of 4.5 cm³, 1.75 cm³ and 0.25 cm³.

Feature Benefits

- Automatic Functions -** Operating simplicity is achieved via PC and a standard web browser or via built-in keypad following prompts or instructions which are automatically presented on the alphanumeric display. Measurements and reporting of results are totally automatic. Continuous self-diagnostics monitor and signal fault conditions that may arise. The transducer is automatically reset to zero prior to each run. Front panel LEDs display the operational status at all times.
- Sample Identification -** In the analysis parameters, alphanumeric may be added to uniquely identify the sample, even from keypad.
- Temperature Readout -** Sample temperature is displayed and printed to ± 0.1 °C. This feature is important for:
 (a) verifying operation at the calibration point or,
 (b) making corrections when analyzing larger quantities of materials whose density varies significantly with temperature.
- Target Pressure -** The measurement of pressure-deformable cellular foams (insulation) is made possible by this feature. A user can conveniently reduce the target pressure typically from around 18 psi to 2 psi. Special versions, the **ULTRAFOAM** and **PENTAFOAM** use optimum low pressure and perform calculations of open and closed cell content with or without ASTM correction for cut cells.
- Sample Purge -** Before analysis, samples are automatically conditioned to remove contaminants and trapped air. The user has a choice to purge by a *continuous flow* or, by a pulse mode suitable for powders. The **PENTAPYC 5200e** purges all five stations simultaneously, thus requiring no more time than for a single analysis. The **ULTRAPYC 1200e** has the added capability to purge the sample by vacuum for a “user selectable” time, (vacuum pump supplied separately).
- Repeat Run Mode -** This feature eliminates the need to reenter the same set of analysis parameters prior to each run. It allows one to quickly change sample weight and sample ID, or simply to rerun a sample by a double keystroke. Greater simplicity and ease of operator is afforded when connected to a PC.
- Useful Statistics -** For three or more averaged measurements the % coefficient of variance and the standard deviation of the volume and the density are reported. This allows a relative comparison of samples from run-to-run and a precise assessment of the absolute variation of the sample being measured.
- NIST Traceability -** Pycnometer volume calibration spheres can be provided with a formal Report of Calibration from the National Institute of Standards and Technology.
- PC Connectivity -** Both Pycnometers feature USB connections to a PC and communication via a standard web browser. Analysis can be started and monitored from the PC, even remotely on a network. Reports can be accessed for analysis on the PC, or can be stored on an external USB flash / thumb drive.
- Balance Port -** Permits interface with configurable analytical balance for automatic transfer of sample weight. Eliminates risk of transcription error.

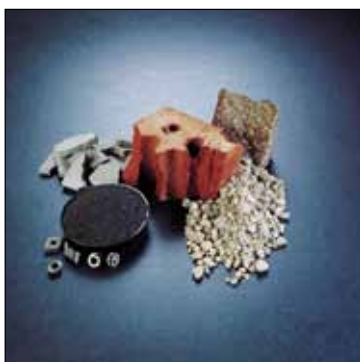

Powder Metals and Construction



Foams and Fibers



Gas Input and Vacuum Fittings



Quantachrome Instruments' corporate headquarters in Boynton Beach, Florida.

Quantachrome®

Renowned innovator of ideas for today's porous materials community.

For 45 years, Quantachrome's scientists and engineers have revolutionized measurement techniques and designed instrumentation to enable the accurate, precise, and reliable characterization of powdered and porous materials:

- Adsorption/Desorption Isotherms
- Surface Area Measurement
- Pore Size Distribution
- Chemisorption Studies
- Water Sorption Behavior
- Mercury Porosimetry
- True Solid Density
- Tapped Density

Not only are Quantachrome products the instruments of choice in academia, but the technology conceived and developed by our expert staff is applied in industrial laboratories worldwide, where research and engineering of new and improved porous materials is ongoing. Manufacturers also rely on porous materials characterization technology to more precisely specify bulk materials, to control quality, and to isolate the source of production problems with greater efficiency.

Quantachrome is also recognized as an excellent resource for authoritative analysis of your samples in our fully equipped, state-of-the-art powder characterization laboratory.



Quantachrome Instruments Application Laboratory.

CORPORATE HEADQUARTERS

Quantachrome Instruments

1900 Corporate Drive
Boynton Beach, FL 33426 USA
Phone: +1 (561) 731-4999
Fax: +1 (561) 732-9888
E-mail: qc.sales@quantachrome.com
www.quantachrome.com

CHINA

Quantachrome Representative Office

M806. Jingbao Garden
183 Andingmenwai Street
Beijing 100011, China
Tel/Fax: +86 10 64400892
www.quantachrome-china.com

EUROPE

Quantachrome UK Limited

Pale Lane Farm, Pale Lane
Hartley Wintney
Hook RG27 8BA, UK
Phone: +44 (0) 125 281 9719
Fax: +44 (0) 125 281 9901
www.quantachrome.co.uk

EUROPE

Quantachrome GmbH & Co. KG

Rudolf-Diesel Str. 12
85235 Odelzhausen, Germany
Phone: +49 (0) 8134 93240
Fax: +49 (0) 8134 932425
www.quantachrome.de

JAPAN

Quantachrome Instruments Japan G.K.

KSP W311 Sakado 3-2-1, Takatsu-ku
Kawasaki-shi, Kanagawa 213-0012
Japan
Phone: +81 44 829 1525
Fax: +81 44 829 1527
www.quantachrome.co.jp



Quantachrome Instruments' quality management system is certified to be in accordance with ISO9001:2008.

WORLDWIDE SALES AND SERVICE

Argentina	Mexico
Australia	Middle East
Austria	Morocco
Bahrain	Netherlands
Belarus	New Zealand
Belgium	Norway
Brazil	Oman
Bulgaria	Pakistan
Canada	Peru
Central America	Philippines
Chile	Poland
China	Portugal
Colombia	Puerto Rico
Croatia	Romania
Cyprus	Russia
Czech Republic	Saudi Arabia
Denmark	Singapore
Egypt	Slovak Republic
Estonia	Slovenia
Finland	South Africa
France	South Korea
Germany	Spain
Greece	Sri Lanka
Hungary	Sweden
India	Switzerland
Indonesia	Taiwan
Ireland	Tanzania
Israel	Thailand
Italy	Turkey
Jamaica	Ukraine
Japan	United Arab Emirates
Jordan	United Kingdom
Kuwait	Uruguay
Latin America	Uzbekistan
Latvia	Venezuela
Lithuania	Vietnam
Malaysia	

Quantachrome
INSTRUMENTS

Serving Porous
Materials and Powder
Characterization
Needs Since 1968



www.quantachrome.com

Trademarks and registered trademarks are the property of their respective owners.