


## Salt-bridge junctions



| Type             | AL100                                | AL110   | AL120                                    |
|------------------|--------------------------------------|---------|--|
| Part no.         | B40A520                              | B40A550 | B40A610                                  |
| Lower diameter   | 12 mm                                | 8 mm    | 8 mm                                     |
| Junction         | Ceramic                              | Ceramic | Ceramic                                  |
| Special features | for Ø 7.5 or 8 mm<br>ref. electrodes |         | for ref. electrodes<br>with ground joint |

During measurement, the salt-bridge junction (also called liquid junction protection tube) prevents the reference electrode from being contaminated by the media and vice-versa. It consists of a glass tube terminated by a liquid junction. Always ensure the filling solution used corresponds to the application.

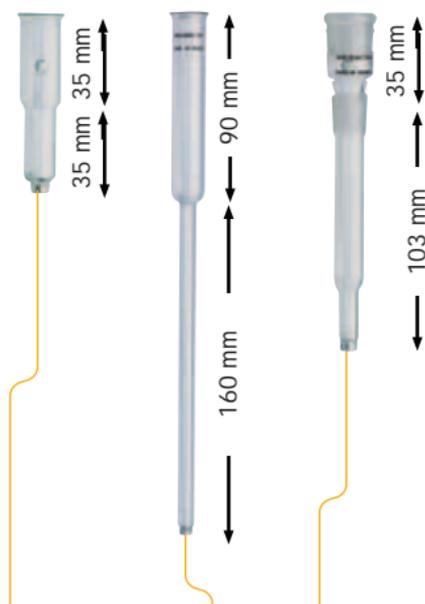

A Hach Company Brand